


TOUCHBeauty

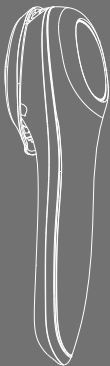

CREAM BOOSTER

SKU24545

INSTRUCTION MANUAL

Safety Instructions:

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

WARNING

- This device is not a toy, keep out of reach of children.
- Please verify the device before you start using it. Do not use if there is any damage.
- Keep the USB cable and the USB charging port dry.
- Do not use the device while charging.

CAUTION


- The following persons should not use this device: pregnant women or while nursing or during menstruation; persons allergic to light; persons with facial nerve disorders; persons who have fractures or implants with metal, silicone, or plastic; persons with open wounds, chapped, or swollen skin on the area of use; persons with implanted electronic medical devices or metal prosthesis; persons with other severe skin diseases.
- This device is not intended for use by persons (including children) with reduced physical, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the device by a person responsible for their safety.
- Immediately stop using the device and consult a doctor if you feel any discomfort.
- Do not look directly at the light source to avoid eye discomfort.
- Keep the device out of reach of children.
- Regularly clean the device to extend its lifetime.
- Do not share the same device with others.
- Do not attempt to disassemble, repair, or modify the device by yourself to avoid accidents.
- This device has a lithium battery. Do not place or store the device for a long period under a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

Device Specifications:

- Name: Cream Booster
- SKU24545
- Model No. TB-1385
- Input: 5V 0.5A
- Type of Battery: Li-ion 3.7V 500mAh
- Executive Standard: Q/XM 1002-2020 Q/XM 1021-2020

Device Description:

- ① Yellow Light
- ② Massage Head
- ③ Cotton Pad Fixing Ring
- ④ Light/Charging Indicator
- ⑤ Ion- Indicator
- ⑥ Ion+ Indicator
- ⑦ Yellow Light Button
- ⑧ Ion- Button
- ⑨ Ion+ Button
- ⑩ Silver Touch Bar
- ⑪ Charging Port
- ⑫ USB Cable

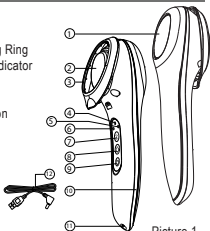

Picture 1

Preparation Instructions:

※ This device battery is low before use. Please charge the device with the USB cable before you start using it.

Charging the Device

- Make sure the device is turned off.
- Connect the USB cable to the device's USB charging port located on the bottom of the device. Then plug the adapter into a wall socket (adapter is not included). (Picture 2)
- While charging, the charging indicator lights up red. When the battery is fully charged, the charging indicator lights up blue. Charging the device takes approximately 3 hours.
- When the battery is fully charged, the device can be used consecutively for about 120 minutes.

Caution: Do not attempt to replace the device's battery.



Picture 2

How to Use:

To ensure the correct operation of this device, please read the instructions carefully before use.

1. Button Instructions

❖❖❖ Light Button

Press once: to turn on the device and the Yellow Light Mode, the Light Indicator lights up red continuously.

Press a second time: to change to Intermittent Yellow Light Mode, the Light Indicator lights up red intermittently.

Press a third time: to turn off the device.

❖❖❖ Ion- Button

Press once: to turn on the device and the Negative Electrode Mode, the Ion- Indicator lights up blue continuously.

Press a second time: to change the function to Intermittent Negative Electrode Mode, the Ion- Indicator lights up blue intermittently.

Press a third time: to turn off the device.


❖❖❖ Ion+ Button

Press once: to turn on the device and the Positive Electrode Mode, the Ion+ Indicator lights up green continuously.

Press a second time: to change the function to Intermittent Positive Electrode Mode, the Ion+ Indicator lights up green intermittently.

Press a third time: to turn off the device.

2. Recommended Use

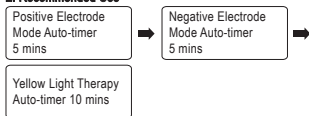

※ It is recommended to use the intermittent mode for people with sensitive skin or for first-time users. Each mode has an auto-timer that turns the device off automatically after a specific period of time.

3. Using Steps

※ Remove makeup and cleanse the skin before using the device.

1) Open the Cotton Pad Fixing Ring, place the cotton pad on the massage head and close the Cotton Pad Fixing Ring to fix it. Apply toner onto the cotton pad. With the device in Positive Electrode Mode, ensure your fingers are in contact with both sides of the Silver Touch Bars. Place the Massage Head on your face and move it gently in a circular motion from the center to the sides of the face. After you finish this mode, remove the cotton pad. (Picture 3)



Picture 3

2) Change the function to Negative Electrode Mode. Apply moisturizing cream/serum to your face and ensure your fingers are in contact with both sides of the Silver Touch Bars. Place the Massage Head on your skin and move it gently in a circular motion from the center to the sides of the face to increase skincare products absorption. (Picture 4)



Picture 4

3) After you finish the Negative Electrode Mode, change the function to Yellow Light Therapy. Keep the device 2 to 3 cm away from your skin and evenly irradiate each area of the face. (Picture 5)



Picture 5

4) After use, turn off the device and clean it.

Cleaning & Maintenance:

※ To ensure the best performance of this device, clean it with a clean slightly wet cloth before and after each use.

1. Make sure the device is power off before you start cleaning.
2. Do not immerse the device in water.
3. Do not use steel wool, abrasive cleansers, chemicals or corrosive liquids to clean the device.
4. Do not use lubricants, detergents, polishing agents or air fresheners on this device.

Storage:

This device has a lithium battery. Do not place or store the device for long periods under a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

Recycling:

The built-in battery contains substances that may pollute the environment. When disposing of this device, do not put it in the general household waste. Please dispose of this device in an official recycling point to contribute to environmental protection.

Warranty & Support:

1. Troubleshooting

Problem	Possible Cause	Solution
The device does not work.	The battery is running low or out of power.	Charge the device as guided.
Other problems.		Contact your local dealer.
Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local dealer.		

2. After-Sales Service

Please visit our website www.touchbeauty.com or contact our customer service team customer@touchbeauty.com if you need to learn about our products or if you need any information. You can also contact an affiliate dealer in your country.

MANUFACTURER: TOUCHBEAUTY BEAUTY & HEALTH
(SHENZHEN) CO., LTD

ADDRESS: 29/F, Block B, Tiley Central Plaza, No.195 Haide
3rd Road, Nanshan District, Shenzhen, China