

TOUCHBeauty

SKU24550

Eye Beauty Device
Instruction Manual



- ※ Important Notice: This manual is available in two languages: English and Simplified Chinese. If any dispute or confusion arises caused by the translation, the English version shall prevail.
- ※ To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

WARNING

- ⊘ This device is not a toy. Keep it out of reach of children.
- ⊘ Please verify the device before you start using it. Do not use it if there is any damage.
- ⊘ Keep the charging cable and charging port dry.
- ⊘ Do not use the device while charging.

CAUTION

- ⊘ The following persons should not use this device: persons with allergies; persons with facial nerve disorders; persons with severe skin diseases on the face; etc.
- ⊘ Do not use the device while wearing contact lenses.
- ⊘ This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the device by a person responsible for their safety.

- ⊗ Immediately discontinue use and consult a doctor if you feel any discomfort.
- ⊗ Keep the device out of reach of children.
- ⊗ Do not share the same device with others.
- ⊗ Regularly clean the device and the accessories to extend its lifetime.
- ⊗ Do not attempt to disassemble, repair, or modify the device by yourself to avoid accidents.
- ⊗ This device has a lithium battery. Do not place or store the device for a long time in humid or hot places. If you will not be using the device for an extended period, charge the device before storing it. Do not throw the device into a fire or put it in the general household waste. Please dispose of this device at an official recycling point.

Device Specifications


TB-1883

- ⊗ Name: Eye Beauty Device
- ⊗ Model No.: TB-1883
- ⊗ Input: 5V==0.5A
- ⊗ Type of Battery: Li-poly 3.7V 800mAh
- ⊗ Executive Standard: Q/XM 1002-2020
Q/XM 1015-2020

Device Description

TB-1883

- ① Eye Massage Head
- ② MASSAGE Mode Indicator
- ③ LIFT Mode Indicator
- ④ Intensity Level Indicator Light
- ⑤ Charging Indicator
- ⑥ Power/Mode Button
- ⑦ Type-C Charging Port
- ⑧ Charging Cable

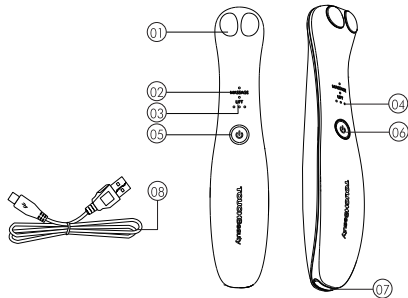

Picture 1

※ This device's battery is low before use. Please charge the device with the type-c charging cable before using it for the first time.


※ The device has a low battery reminder function. When the battery is low, the charging indicator light flashes 10 times with 5 low-battery sounds. Please charge the device before it runs out of battery.

🔌 Charging the Device:

- Make sure the device is turned off.
- Connect the charging cable to the device's charging port located on the bottom of the device. Then plug the adapter into a wall socket (adapter not included). (Picture 2)
- When charging, the charging indicator light flashes. When fully charged, the charging indicator light is always on. Charging time is approximately 3 hours.
- A full battery allows for approximately 60 minutes of continuous use. You can use it 7-8 times, 8 minutes each time.

Caution: Do not attempt to replace the device's battery.

Picture 2

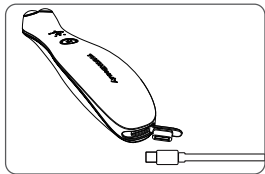

To ensure the correct operation of this device, please read the instructions carefully before use.

Device Operation

1) Mode Instructions

- The **MASSAGE** mode combines 42°C hot compress with constant temperature + sonic vibration to help to soothe the skin.

You can choose to use this mode with or without a skincare product. If you use this mode with a skincare product, choose one that fits your skin needs.

- The **LIFT** mode combines EMS microcurrent + 42°C hot compress with constant temperature to help to lift and tighten the skin, smooth the appearance of fine lines around the eyes.


Use this mode with eye essence, eye cream or gel to boost efficiency and conduct current. We do not recommend using this mode on dry skin.

2) Button Instructions

① **Turn-ON/OFF:** Press and hold the power/mode button to turn ON/OFF the device.

② **Mode Switch:** Short press the power/mode button to switch between modes and intensity levels. The intensity level indicator on the device is "MASSAGE level 1" - "MASSAGE level 2" - "MASSAGE

level 3" - "LIFT level 1" - "LIFT level 2" - "LIFT level 3". Switch between 6 modes in a cycle. The MASSAGE mode has 3 intensity levels that control the vibration intensity, while the LIFT mode has three levels that control the microcurrent intensity. You can adjust the intensity level from low to high. (Picture 3)

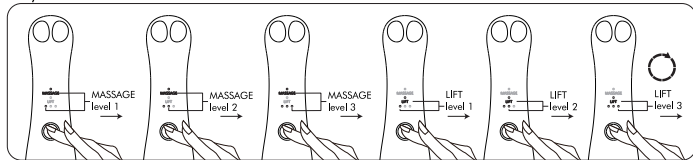

Picture 3

③Quick Mode Switch: When you turn ON the device, it automatically starts the MASSAGE mode by default. Press the power/mode button twice in quick succession to switch between the modes directly.

※Time Reminder

If there is no button operation with the power on, the device automatically turns OFF after 6 minutes.

👁Using Instructions

※It is not recommended to use the device on the eyelids to prevent hurting the eyes. Do not use this device if your eyes are inflamed or swollen.

※Eye Massage

- Press and hold the power/mode button to turn on the device and the MASSAGE mode. Select an intensity level according to your skin needs and preference. Gently press the massage head onto the skin, slowly glide it or apply light pressure around the eye contour. We recommend using the MASSAGE mode for 1 min on each eye.(Picture 4)



Picture 4

※Eye Wrinkles

- Before you start using the LIFT mode, use the MASSAGE mode to help to wake up the skin around the eyes and boost the absorption of skincare products. Apply a layer of eye cream/essence to the skin around the eyes.

Switch to LIFT mode and select the microcurrent intensity level according to your preference and skin needs. Gently press the massage head onto the skin and slowly glide it from the inner corner of the eye towards the temple area. Afterward, glide the massage head along the outer corner of the eye in repeating lifting strokes. We recommend using the LIFT mode for 1 min on each eye. (Picture 5)

- Use daily for better results.



Picture 5

※ After using the device, turn it OFF and follow your skincare routine.

Cleaning & Maintenance

TB-1883

- ✕ Before and after every use, please clean the device body to ensure the best performance of this device.
- 🚫 Make sure the device is turned off before you start cleaning.
- 🚫 Clean the device with a clean, slightly wet cloth. Do not immerse the device in water.
- 🚫 Do not use alcohol, steel wool, abrasive cleansers, chemicals, or corrosive liquids to clean the device.
- 🚫 Do not use lubricants, detergents, polishing agents, or air fresheners on this device.

Storage

TB-1883

This device has a lithium battery. Do not place or store the device for a long time in humid or hot places. If you will not be using the device for an extended period, charge the device before storing it.

Recycling

TB-1883

The built-in battery contains substances that may pollute the environment. When disposing of this device, do not put it in the general household waste. Please dispose of this device in an official recycling point to contribute to environmental protection.

Troubleshooting

Problem	Possible Cause	Solution
The device does not work.	The battery is running low or out of power.	Charge the device as guided.
The current sensation feeling is not strong.	Low conductivity.	Apply conductive essence/gel.
Other problems.	Requires specific analysis.	Consult your local dealer.
Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local dealer.		

☹ After-Sales Service

Please visit our website www.touchbeauty.com or contact our customer service team customer@touchbeauty.com if you need to learn about our products or if you need any information. You can also contact an affiliate dealer in your country.

☹ MANUFACTURER: TOUCHBEAUTY BEAUTY & HEALTH (SHENZHEN) CO., LTD

☹ ADDRESS: 29/F, Block B, Tiley Central Plaza, No.195 Haide 3rd Road,
Nanshan District, Shenzhen, China