├ TOUCHBeauty

SKU24555

Ultrasonic Beauty Device Instruction Manual



Safety Instructions TB-1887

Important Notice: This manual is available in two languages: English and Simplified Chinese. If any dispute or confusion arises caused by the translation, the English version shall prevail. To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

AWARNING

- This device is not a toy, keep out of reach of children.
 Please verify the device before you start using it. Do not use if there is any damage.
- 🚭 If used in a dry environment, handling this device may generate static electricity. This is a normal phenomenon. It does not affect the use of the device.
- Keep the USB cable and the charging base dry. Do not use this device on dry skin.
- Do not use the device repeatedly in the same area over a short time.
- Do not use this device with strong exfoliating products that are abrasive or contain granules.
 Do not use the device while charging.

ACAUTION

The following persons should not use this device: pregnant women or while nursing or during menstruation; persons with allergies; persons with facial nerve disorders; persons with cardiac

dysfunctions, acute disease, or malignant tumors; persons who have fractures or implants with metal, silicone, or plastic; persons wearing electronic medical equipment; persons who have taken steroid hormones medication for a long time; persons with open wounds, chapped or swollen skin

in the area of use; persons with severe skin diseases on the face. a Do not use this device in the following areas: eyes and ears; sensitive areas with wounds, eczema, skin inflammation, etc.; areas where surgery has been performed; areas with implanted silicone, metal, resin, etc.; areas with acne, skin ulcers or other skin diseases. This device is not intended for use by persons (including children) with reduced physical, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the device by a person responsible for their safety.

Immediately stop using the device and consult a doctor if you feel any discomfort. Seep the device out of reach of children.

Regularly clean the device to extend its lifetime.

Do not share this device with others.

Do not attempt to disassemble, repair, or modify the device by yourself to avoid accidents. Tie your hair before using the device and remove jewelry, to prevent injury during use.

This device has a lithium battery. Do not place or store the device for long periods under a hot or

humid environment. If you do not intend to use the device for a long time, store it after it is fully charaed.

Device Specifications


TB-1887

- Name: Ultrasonic Beauty Device
- Model No.: TB-1887 Input: 5V==−/0.5A
- Power: 3.7V === /2.4W
- Type of Battery: Li-ion 3.7V 800mAh
 Executive Standard: Q/XM1024-2020

Device Description

TB-1887

- Silicone Protective Cap (II) Cleaning Cloth
- @Ultrasonic Spatula Blade @ Intensity Level Button
- Power/Mode Button
- Sensor Strip
- Device Contact Points
- O Charging Base Contact Points
- ® Base Charging Port
- @ USB Cable

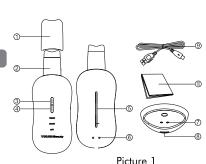



Preparing the Device for Use TB-1887

*This device battery is low before the first use. Please charge the device with the USB cable before you start using it.


The device has a low battery alarm system. When the battery is low, the intensity indicators will flash three times, and the device will automatically turn off.


Charging the Device

- · Make sure the device is turned off.
- Connect the USB cable to the wireless charging base and to an adapter (adapter not included).
- Place the device on the charging base, and plug the power adapter into the wall socket. (Picture 2)

 While charging, the intensity indicators light up in a cycle. When the battery is full, the intensity indicators light up continuously. Charging the device takes approximately 4 hours.

Note: Do not attempt to replace the device's battery.

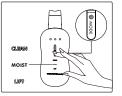

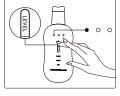

Picture 2

How to Use

TB-1887

To ensure the correct operation of this device, please read the instructions carefully before use.

Device Operation

Power/Mode Button

Press and hold the power/mode button to turn on the device and the "CLEAN" mode. Short press the power/mode button to change between modes "CLEAN" —> "MOIST" —> "LIFT". When in LIFT mode, press the power/mode button to turn off the device. To turn off the device in any mode, press and hold the power button. (Picture 3)

*This device has an auto-timer. It will automatically shut down after 5 minutes of use.

Intensity Level Button

After turning on the device, the intensity indicator lights up the first intensity level. Short press the intensity level button to change between intensity levels in a cycle "intensity 1", "intensity 2", "intensity 3". (Picture 4)






Recommended Use


deep cleansing the skin. (Picture 5)

While using the device, hold the body so your fingers/hand touch the sensor strip on the back of the device.

the device.

CLEAN Mode:
 We recommend using the 'CLEAN' mode 1-2 times a week. Do not use this mode for more than 5 minutes each time. Remove makeup and cleanse the face before you start using the device. Dampen the skin with water or toner. Turn on the device and select 'CLEAN' mode. Adjust the intensity level according to your needs. Place the front (angled) side of the ultrasonic spatula blade on the skin at a 30° angle. Slowly move the spatula along facial contours as illustrated. Using this mode helps to

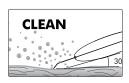

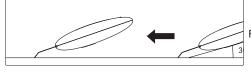


Picture 5

MOIST Mode:

We recommend using the 'MOIST' mode 3-5 times a week. Do not use this mode for more than 5 minutes each time. Apply your daily skincare products. Select the 'MOIST' mode and adjust the intensity level. Place the back (flat) side of the ultrasonic spatula blade on the skin at a 30° angle. Slowly move the spatula along facial contours as illustrated. Using this mode helps to boost the absorption of skincare products. (Picture 6)



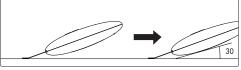


Picture 6

LIFT Mode:

We recommend using the 'LIFT' mode 3-5 times a week. Do not use this mode for more than 5 minutes each time. After using the 'MOIST' mode, select the 'LIFT' mode and adjust the intensity level. Place the backside of the ultrasonic spatula blade on the skin at a 30° angle. Slowly move the spatula along facial contours as illustrated. Using this mode helps to lift and firm the skin. (Picture 7)




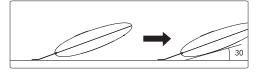

Picture 7

Areas of Use

• Use the device in the following areas (refer to below picture):

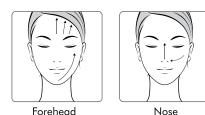





※Turn off the device after use and clean the ultrasonic spatula blade with the cleaning cloth. Afterward, follow up on your daily skincare routine.

Cleaning & Maintenance

TB-1887

To ensure the best performance of the device, clean it with a clean slightly wet cloth before and after each use. Make sure the device is power off before you start cleaning.

Do not immerse the device in water. Use a clean cloth or a slightly wet cloth to clean the device body and the ultrasonic spatula blade.

Use a dry cloth to clean the charging base. If some areas are difficult to clean, wipe it with a slightly wet cloth.

Do not use lubricants, detergents, polishing agents, or air fresheners on this device.

Storage

TB-1887

This device has a lithium battery. Do not place or store the device for long periods under a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

Recyclina

TB-1887

The built-in battery contains substances that may pollute the environment. When disposing of this device, do not put it in the general household waste. Please dispose of this device in an official recycling point to contribute to environmental protection.

Warranty & Support

TB-1887

Troubleshooting

| Problem | Possible Cause | Solution |
|--|---|---|
| Press the power/mode button, the device does not work. | The battery is low our out of power. | Please charge the device according to the instructions. |
| Unable to produce water mist. | The battery is low. | Please charge the device according to the instructions. |
| Unable to charge the device. | The contact points in the device are not in contact with the charging base correctly. | Adjust the device position until the contact points are in contact. |
| Other problems. | | Consult your local dealer. |

Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local dealer.

After-Sales Service Please visit our website www.touchbeauty.com or contact our customer service team customer@touchbeauty.com if you need to learn about our products or if you need any information. You can also contact an affiliate dealer in your country.

■ MANUFACTURER: TOUCHBEAUTY BEAUTY & HEALTH (SHENZHEN) CO., LTD.

ADDRESS: 29/F, Block B, Tiley Central Plaza, No.195 Haide 3rd Road, Nanshan District, Shenzhen, China