
# **TOUCH**Beauty

SKU24559

Multi-Therapy Device Instruction Manual

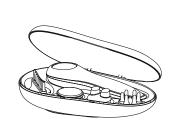

## Safety Instructions TB-1713

※Important Notice: This manual is available in two languages: English and Simplified Chinese. If any dispute or confusion arises caused by the translation, the English version shall prevail.

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

#### **<b>∆**WARNING

The following persons should not use this device: pregnant women or while nursing; children; persons who have a pacemaker or other electronic medical equipment; persons with allergic reactions; persons with skin diseases (skin blemishes, burns, etc.); persons who epilation or plastic surgery was performed on the areas of use; on the eyelids and throat; while wearing makeup.

1

- This device has different heads suitable to use on the face (forehead, cheeks, and chin) and scalp. Do not use it on the evelids. If you are light sensitive, prone to rash or allergic reactions, do not use the light therapy mode.
- Do not use the light therapy mode if the treatment area has infections, eczema, burns, inflammation of hair follicles, open lacerations, herpes, wounds or lesions, and hematomas.
- Persons with diabetes, lupus erythematodes, porphyria or congestive heart disease, bleeding disorder, history of immunosuppressive disease (including HIV infection and
  - AIDS) should not use this device.
  - If your skin is currently being treated or has recently been treated with Alpha-Hydroxy

Do not use this device if you are undergoing radiotherapy or chemotherapy. Do not use this device if you take any of the following medication:

● If you have taken any form of Accutane A or Roaccutane A in the past six months. This

Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.

treatment makes skin susceptible to laceration, trauma, and irritation.

- If you are taking photosensitizing agents or any medication, check the patient information leaflet before you start using this device. Do not use this device if the information leaflet specifies that it can cause photo-allergic and photo-toxic reactions, or avoid exposure to sunlight while taking the medicine.
- If you take anticoagulant medication, including large doses of aspirin, only use the device one week after you stop the treatment.
- Do not repair, disassemble (except interchangeable heads) or modify the device, as improper handling may damage it. Do not use this device for any other purpose.
- Do not look directly at the light when using the light therapy head as it may cause discomfort. Use protective goggles to protect the eyes.
- This device has a lithium battery. Do not place or store the device for a long period under a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged. Do not throw the device into a fire or put it in the general household waste. Please dispose of this device in an official recycling point.

## Device Specifications TB-1713

- Name: Multi-Therapy Device
- Model No.: TB-1713
- Input: 5V == /0.5A
- Power: 3.7V === /1W
- Type of Battery: Li-ion 3.7V 500mAh
- Executive Standard No: Q/XM1002-2020 Q/XM1021-2020

# **Device Description**


TB-1713

① Device Body

**4 USB Cable** 

② Connection Point③ Power Button

- (§) Cleaning Cloth (§) Storage Box
- 4 ION+ / / Microcurrent Button
- (A) ICIN+ / / MICIOCUITEIII BUILI
- ⑤ LED Button
- Conductive Bars
- 7 Charging Port
- ® Light Therapy Head
- ① Cotton Ring
- 11 Microcurrent Facial Head
- Scalp Massage Head
- 13 Microcurrent Scalp Head

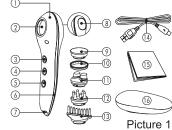

# Certification Marks TB-1713



## Preparation Instructions | TB-1713

- Charging the Device
- % Charge the device according to its charging system. Charging the device takes approximately 2.5 hours. When the battery is fully charged, the device can be used consecutively for about 1.5 hours.
- The light indicator flashes red when charging and it turns blue when the device battery is full.
- $\divideontimes$  If the battery is low, the microcurrent indicator flashes 3 times and the device automatically shuts down.

★ The device does not work while it's charging. Connect the charging cable as illustrated. (Picture 2)



Picture 2

- Disassemble/Assemble the Heads
- Disassemble: Hold the device body with one hand, with the other hand rotate the accessory head to the left and gently pull it to remove. (Picture 3)





Picture 3

 Assemble: Hold the device body with one hand, connect the accessory head with the device body, and rotate to the right until it locks the position. (Picture 4)

#### Button Instructions Power Button:

\*You can turn on the device even if the head is not installed.

Press the power button, after a beeping sound the indicator lights up white, and the device enters the standby mode. If there is no operation, the device automatically turns off after 30 seconds.

#### Ion+/- / Microcurrent Button:

%If the corresponding head is installed correctly, press the ion+/- / microcurrent button to turn on the respective function.

\*The device automatically recognizes the function of the accessory head installed. If you change the accessory head, the ion+/- / microcurrent button will display the corresponding mode.


- After installing the ion+/- beauty head, press the ion+/- / microcurrent button, the indicator lights up green after a sound and the device starts the lon+ (cleansing) mode. Press the ion+/- / microcurrent button again, the indicator lights up yellow after a sound and the device changes to lon- (moisturizing) mode. The two modes change in a cycle.
- After installing the microcurrent scalp head or the microcurrent facial head, press the ion+/-/microcurrent button. After a sound, the indicator lights up green, and the device starts the microcurrent low-intensity mode. Press again the lon+/-/microcurrent button, after a sound the indicator lights up yellow, and the device switches to microcurrent high-intensity mode.

#### LED Button:

- The LED light therapy head is incorporated into the device. Press the LED button to start the light therapy mode.
- Turn on the device and short press the LED button. After a sound, the indicator lights up red
  and the red-light mode starts. Press the LED button again, after a sound, the indicator lights
  up yellow and the device changes to flashing yellow light therapy mode.

## Using Time

Facial Care:

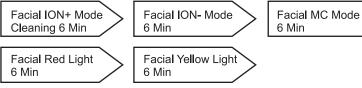

Head Care:

Head ION+ Mode Cleaning 6 Min Head ION- Mode 6 Min Head MC Mode 6 Min

This device is suitable for daily use.

## How to Use TB-1713

To ensure the correct operation of this device, please read the instructions carefully before use.

- Daily Facial Care
- ※ Use the device on clean skin, be sure to remove makeup and impurities from the skin before use.
- $\divideontimes$  Do not look directly at the light when using the light therapy head as it may cause discomfort on the eyes.

#### Ion + Mode

 Install the facial beauty head, place a cotton pad on the head surface and press the cotton ring down over the cotton pad. Soak toner into the cotton pad. (Picture 5)



- Press the power button and the ion+/- / microcurrent button in sequence to start the cleansing mode. This mode has sonic vibration.
- Hold the device with both sides of your hand, so it touches the conductive bars. With the ion+/- beauty head touching the skin, gently move it in a circular motion from the center to the sides of the face. (Picture 6)



Picture 6

#### Ion - Mode

 After using the ion+ cleansing mode, remove the cotton pad and cleanse the face with warm water. Switch to Ion- moisturizing mode. • Apply lotion/serum to the face. Hold the device with both sides of your hand. Place the surface of the accessory head on your face and gently move it in a circular motion from the center to the sides of the face to promote the absorption of skincare products. (Picture 7)



Picture 7

#### Facial MC Mode

- After using the ION- mode, install the microcurrent facial head.
- Press the power button and the ion+/- /microcurrent button in sequence. The device starts
  the microcurrent low-intensity mode and sonic vibration. Press again the ion+/- /MC button
  to start the microcurrent high-intensity mode.
- Slide the surface of the accessory head from the lower jaw to the cheekbone or from the corner of the mouth to the ear. Select the intensity of the microcurrent mode, after a

short sound, move the surface of the accessory head as illustrated in picture 8. Repeat the procedure 2 times.





Picture 8

• Around the mouth: Select the microcurrent intensity, after a sound, place the surface of the accessory head on the skin and move it as illustrated in picture 9. Repeat the procedure 2 times.

Picture 9

• Cheeks: To use the accessory head on the cheeks, move it from the corner of the mouth to the ear. After a sound, slightly move the massage head as illustrated in picture 10. Repeat the procedure 2 times. (Picture 10)

Picture 10

• Around the Eyes: Massage as illustrated by the arrows in picture 11 and repeat the procedure 2 times.

Picture 1

• Forehead: Close the eyes, move the surface of the accessory head slowly along the forehead up to the hairline as illustrated in picture 12, repeat the procedure 2 times.



Picture 12

### Red/Yellow Light Therapy

\*\* Test if your skin is light sensitive prior to use: Use the light therapy mode on the inside of your forearm for 3 minutes. If after 3 minutes of use the area has turned red and the redness lasts for more than 2 hours, your skin is light sensitive. Please stop using the light therapy mode.

☼ Do not look directly at the light therapy source to avoid discomfort on the eyes.

• Press the power button and the LED button in a sequence to turn on the red-light mode. Press the LED button again to change to yellow light therapy mode (intermittent mode).

• Place the light head within a 2-3cm distance from the skin to evenly irradiate the surface of the skin. The light therapy head can also be in contact with the skin while using this mode. Use the device for 2 minutes on the left side of the face, 2 minutes on the right side of the face, and 2 minutes on the forehead. Every 6s the device emits a sound reminding you to change the treatment area. (Picture 13)

## Daily Head Care lon + Mode

- Install the scalp massage head.
- Press the power button and the ion+/-/microcurrent button to start the ion+ cleaning mode. This mode includes sonic vibration.
- Comb the scalp to clean it as recommended in picture 14 (order 1-3).

Picture 13







Picture 14

#### Ion - Mode

- Press the ion+/-/microcurrent button again to switch to ion- moisturizing mode. This mode includes sonic vibration.
- Apply preferred scalp product on the scalp.
- Comb the scalp to increase the absorption of the product applied, use as recommended in picture 15 (order 1-3).







Picture 15

#### Microcurrent Mode

- Install the microcurrent scalp head, press the power button and the ion+/-/microcurrent button in a sequence to start the microcurrent mode. Every 6s the device emits a sound reminding you to change the treatment area.
- Apply scalp care product.
- Comb the scalp and massage it in the order 1-3 as recommended in picture 16.







Picture 16

## Cleaning & Maintenance TB-1713

- After using the device, clean the accessory heads and the device body with a clean slightly wet cloth. Do not clean directly with water.
- Do not use steel wool, abrasive cleansers or corrosive liquids to clean the device.

## Storage TB-1713

This device has a lithium battery. Do not place or store the device for long periods under a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

# Recycling TB-1713

The built-in battery contains substances that may pollute the environment. When disposing of this device, do not put it in the general household waste. Please dispose of this device in an official recycling point to contribute to environmental protection.

# Warranty & Support TB-1713

## Troubleshooting

Problem	Possible Cause	Solution
Device is not working properly.	The accessory head is not correctly installed.	Reinstall the head correctly as instructed.
	Battery is out of power.	Charge the device as guided.
Microcurrent is weak.	Skin is too dry.	Use toner or serum.
Other problems.		Please consult your local dealer.

Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local retailer.

After-Sales Service
Please visit our website www.touchbeauty.com or contact our customer service team customer@touchbeauty.com if you need to learn about our products or if you need any information. You can also contact an affiliate dealer in your country.

MANUFACTURER: TOUCHBEAUTY BEAUTY & HEALTH (SHENZHEN) CO., LTD

ADDRESS: 29/F, Block B, Tiley Central Plaza, No.195 Haide 3rd Rood, Nanshan District, Shenzhen, China